

LAFANTA ELECTRO HYDRAULIC OPERATION THEATER TABLE

We offer electro Hydraulic operation theater table powered by European Electro Hydraulic System having capacity to lift 450 kg weight.

Special features

1. All positions operated through cable remote.
Electric Floor locking . Auto leaving (zero position)
2. Built in kidney bridge
3. Three powering options
 - a. 230V AC mains supply
 - b. In build re-chargeable battery.
 - c. Paddle system can be operated manually in event of mains and battery and power failure.
4. Switch can be select three powering options – mains, battery and power off.
5. Hand held remote control and manually foot operates paddle control.
6. Table top slide 300mm enabling excellent performance for c- arm & x –ray.
7. Antistatic PU waterproof & antibacterial matters.
8. Suitable for orthopaedic , urology ,general & Gynaec surgeries .
9. ISO & CE mark.
10. Built in intelligent battery charging facility driven by independent power source.
11. German make sturdy poly urethane castors are used for scratch free movement in modular operation theaters.
12. Sturdy S.S 304 unbreakable base used for preventing water & dust entering inside.



Technical Data:

Top dimension	L 1900*W 520 mm
Height adjustment	720mm-1050mm`
Trendelenburg /Reverse	30°/25°
Lateral tilt	20°/20°
Kidney elevator	130 mm
Back rest (UP/Down)	80°/25°
Leg Rest (UP/Down)	15°/90°
Head Rest (UP/Down)	20°/60°
Power Supply	24V DC
Battery Backup	2-3 Hours (optional)
Patient Weight Capacity	450 kg
Mains Voltage	AC 220 V±10%50 HZ

Electrometric Operation table

Special features

Electrical operated surgical table designed to support Cardiology, Orthopedic, Neurology, Gynecology, Urology and ERCP Procedures. The table is equipped with a large translucent top providing maximum radiological access & c- arm functions. Fully electrically operated table, five sections, radio- translucent table top with eccentric pillars and Die compressed base. Full length C- ARM compatible, having all surgical positions i.e raising / lowering , trendelenburg , reverse Trendelenburg , lateral tilt, flex/ re- flex & chair position . Interchangeable head leg section (manually). All positions operated through cable remote. Supplied complete with standard accessories and mattress. In built Kidney Bridge and NiCD battery backup.



Model: EOT-9401

Technical Data

Top dimension	L 1900*W 520 mm
Height adjustment	720mm-1050mm`
Trendelenburg /Reverse	30°/25°
Lateral tilt	20°/20°
Kidney elevator	130 mm
Back rest (UP/Down)	80°/25°
Leg Rest (UP/Down)	15°/90°
Head Rest (UP/Down)	20°/60°
Power Supply	24V DC
Battery Backup	2-3 Hours (optional)
Patient Weight Capacity	150 kg
Mains Voltage	AC 220 V±10%50 HZ

C-Arm Compatible Hydraulic Operating Table

Hydraulically operated table , five section , radio- translucent table top with eccentric pillars . full length C-Arm compatible having all surgical position i.e. raising / lowering . Trendelenburg , lateral tilt, flex/re-flex & chair position. Interchangeable head and leg section .all position operated through cranks fitted under the head end of the table in built kidney bridge supplied complete with standard accessories and mattress. This is hydraulic surgical table designed especially for othopedic , Neurology , urology and ERCP procedure. Table is constructed of stainless steel and other high quality materials. The table is equipped with a large translucent top, for C-arm use.



Model: CAH-9402

Special features

- Hydraulic and mechanical operated.
- The base is eccentric for trouble free C-ARM movement.
- Leg & head section are interchangeable.
- Eccentric pillars.
- Stainless fitting
- Detachable head, leg and pelvic section.
- Inbuilt kidney bridge.
- Five sectional radio-translucent.
- Base and cover made of attractive, easy to clean stainless steel.
- Sophisticated mechanics provide smooth, stepless articulation of the table top for precise patient positioning.

Top dimension	L 1905*W 553 mm
Height adjustment	730mm-1050mm`
Trendelenburg /Reverse	30°/25°
Lateral tilt	20°/20°
Kidney elevator	150 mm
Back rest (UP/Down)	80°/25°
Leg Rest (UP/Down)	15°/90°
Head Rest (UP/Down)	20°/60°
Patient Weight Capacity	110 kgs.

HYDRAULIC OPERATION TABLE

Hydraulically operated Table, five sections, stainless steel Table top having all surgical positions i.e. Raising / lowering, trendelenburg, reverse Trendelenburg, lateral tilt, flex / re-flex & chair position. Interchangeable head and leg sections. All positions operated through cranks fitted at the side end of the table. Supplied complete with standard SS accessories & matters. This side End

Control OT Table is mobile, hydraulic designed to support virtually all general surgical procedures like: gynecology and urology .This surgical table features can be opting by its mechanism i.e. Lateral tilt,

Trendelenburg / Reverse Trendelenburg, Chair position and Adjustable height functions.



Model: HDT-9403

Special Feature:

- The five section Stainless steel table top is equipped with removable, antistatic mattresses.
- All the positions are realized by a mechanical driver.
- The table is furnished with hard, robust stainless shell which provides the sense of dirt free & sanitized operations.
- Precision made hydraulic system with a single lever action hydraulic life controlled by a foot pedal.
- Inbuilt kidney Bridge and large perinea cut out for drainage.
- Manual floor locking system.

Technical Data

Top dimension	L 1900*W 620 mm
Height adjustment	720mm-1050mm`
Trendelenburg /Reverse	30°/25°
Lateral tilt	20°/20°
Kidney elevator	150 mm
Back rest (UP/Down)	80°/25°
Leg Rest (UP/Down)	15°/90°
Head Rest (UP/Down)	20°/60°
Patient Weight Capacity	110 kgs.

ELECTRO MANUAL OT TABLE

Special features: -

Electrically as well as manually operated table, five sections, stainless steel table top having all surgical positions i.e., raising/lowering, Trendelenburg, reverse Trendelenburg, lateral tilt, flex/re-flex chair position. All positions operated through Hand Remote provided at the side of the table. Supplied complete with standard SS accessories mattress. This Electrically Control OT Table is mobile, designed to support virtually all surgical procedures like: CARM, Gynaecology, Urology and ENT. This Surgical Table features can be opting by Position and Adjustable height functions Electrically.

Salient Feature:

- Electric Control with Battery Backup.
- Longitudinal Top slide 300mm.
- Additional Nurse Panel control.
- Manual over ride control.
- Eccentric Pillars Stainless Fitting.
- Detachable, head, leg & Pelvic section.
- In built Kidney Bridge.
- Five Sectional Radio-Translucent.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning



Model: MMD9405



Nurse Panel Control



Manual over ride control



Hand Remote

Technical Specifications: -

Length of Table	1900MM
Width of Table Top	520 MM
Minimum Height (without Mattress)	Min. 750 MM
Maximum Height (without Mattress)	Max. 1000MM
Trendelenburg	25°
Reverse Trendelenburg	25°
Lateral Tilt	±20°
Back Section	±60°
Head Section	±60°
Leg Section	0°-90°
Top Sliding	12"


OPHTHALMIC HYDRAULIC OT TABLE

Special Features:

This Ophthalmic Table is specially designed for ophthalmic surgery. Three sectional Table is provided with drainage system in the head section. The table has sturdy base, column made up of stainless steel, antibacterial and is easy to clean. Electromotive height with the help of lightweight handheld control unit.

Salient Features:

- ❖ Stainless Fitting
- ❖ Detachable, head & leg section
- ❖ Three Sectional Radio Translucent Top
- ❖ Base & Cover made of attractive, easy to clean Stainless Steel.
- ❖ Sophisticated mechanics provide smooth, step articulation of the table tip for precise patient positioning.

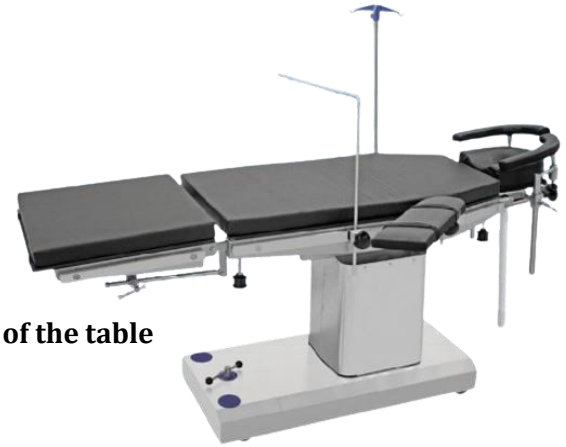

Model No: MMD9406

Technical Data:

Length Of Table	1900mm
Width Of Table Top	520mm
Minimum Height (Without Mattress)	750mm
Leg Section	0°-90°
Battery Backup	Optional



Foot Switch (Optional)

Standard Accessories:

- Anesthetic Screen
- IV Pole
- Hand Rest
- Mattress



Instrument Tray (Optional)


OPHTHALMIC OT TABLE

Special Features:

This Ophthalmic Table is specially designed for ophthalmic surgery. Three sectional Table is provided with drainage system in the head section. The table has sturdy base, column made up of stainless steel, antibacterial and is easy to clean. Electromotive height with the help of lightweight handheld control unit.

Salient Features:

- ❖ Stainless Fitting
- ❖ Detachable, head & leg section
- ❖ Three Sectional Radio Translucent Top
- ❖ Base & Cover made of attractive, easy to clean Stainless Steel.
- ❖ Sophisticated mechanics provide smooth, step articulation of the table tip for precise patient positioning.

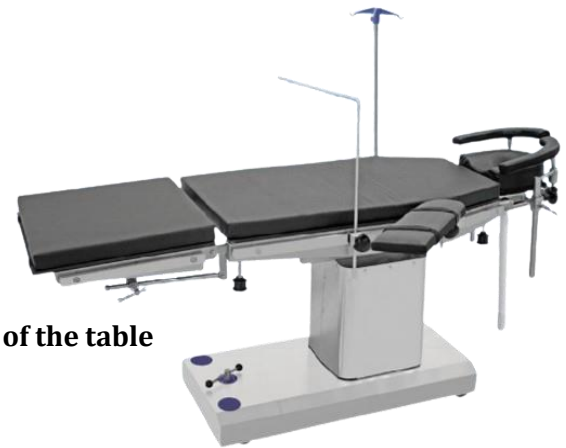

Model No: MMD9406

Technical Data:

Length Of Table	1900mm
Width Of Table Top	520mm
Minimum Height (Without Mattress)	750mm
Leg Section	0°-90°
Battery Backup	Optional



Foot Switch (Optional)

Standard Accessories:

- Anesthetic Screen
- IV Pole
- Hand Rest
- Mattress



Instrument Tray (Optional)